

-80/-40°C

REMEMBER TO GET A QUOTE ON OPTIONS


LOW NOISE LEVEL

ULTRA LOW TEMPERATURE FREEZERS

# MUF 40

MUF (Multi Ultra Freezer)

Whether you need to use it under the counter in your laboratory, or just have some samples stowed always in your private office, the MUF40, with its 35 Liters capacity, is perfect for a quick reach whenever you need.

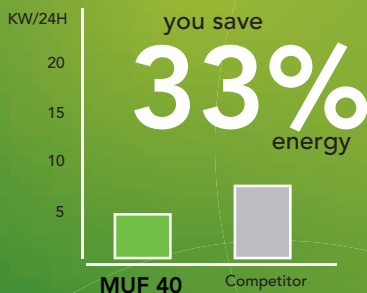

OPTIONAL CASTORS



## MUF 40

- Single compressor
- Low energy consumption
- Low noise level
- Low heat dissipation
- Key lock on door

## Controller features

- Digital display
- Approx. 48 hours battery back up for alarms and temperature display in case of power cut
- Visual and acoustic alarms
- High/low temp. alarm
- Power failure alarm
- Contact for remote alarm

# Portable -80°C freezer

Keep your eyes on what's important.

The DP-80 ensures direct access to your samples at your desk. The Cryo Porter is based on sterling technology which means, no compressor, nor refrigerant. This makes it the most environmentally safe product in the Arctiko range.



REMEMBER THE ALUMINIUM BLOCK

## DP-80

- Aluminium blocks for quick freezing of samples
- Digital controller
- Simple & Easy Operation
- Temp. stability +/- 0,2°C
- Safety system for over current high and low temperature
- Adjustable high/low temperature alarm
- AC (230V) and DC (12V) connection cables

Scan me and see more about MUF 40



Model	DP-80	MUF 40
Item code (DAI no.)	DAI 1420	DAI 1410
Dimensions exterior WxDxH (mm)	230x262x390	495x660x680
Net capacity (L)	-	35
Temperature range (°C)	-80/-40	-80/-40
Power supply (V)	230/12	230/110
Frequency (Hz)	50/60	50/60

## ALUMINIUM BLOCKS:

- (DAI 1421) DPB-0296: 96 wells for PCR-tubes 0,2 ml
- (DAI 1422) DPB-0548: 48 wells for PCR-tubes 0,5 ml
- (DAI 1423) DPB-1548: 48 wells for Sample-tubes 1,5 ml
- (DAI 1424) DPB-1840: 40 wells for Vials, Cryotubes 1,8 ml
- (DAI 1425) DPB-2040: 40 wells, Cryotubes 2,0 ml

We reserve the right to change specifications without notice. Subject to confirmation and availability.