
CO₂ Incubators

High Heat Decontamination SC06AD





The SHEL LAB SCO6AD is designed to stop microbial contamination caused by mycetozoa, yeast, viruses, and bacteria and the range of other microorganisms that thrive in incubator environments.

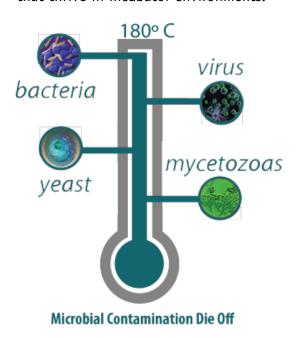

PROTECT YOUR SAMPLES WITH HIGH HEAT DECONTAMINATION!

The model SCO6AD features a dry heat decontamination cycle, that maintains 180°C for 120 minutes. This industry best time and temperature relationship satisfies all global standards for decontamination.

This is decontamination at its fastest easiest and most effective - it is not necessary to remove the IR CO₂ sensor prior to activating the decontamination process and we feature the shortest cycle time on the market. This is a more convenient approach and eliminates potential damage to the sensitive IR sensor.


Other features of the SCO6AD include a USB interface for software communication, precise temperature control microprocessor and an independent over temperature safety controller.

CONTAMINATION CONTROL

Extensive use of copper in the CO2 sample port, heated CO2 feed line and humidity pan with a copper SL decontamination token, adds reassurance that foreign microbes will not affect test results. Cleanup is a breeze with the all stainless steel chamber, and autoclavable door gasket. Optional copper shelves are available for even more contamination control.

Infrared (IR) CO2 Sensors

For the fastest CO₂ recovery and most stable performance, the SCO6AD features IR sensors.

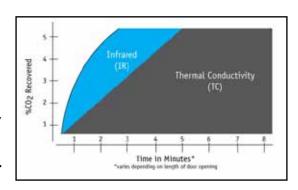



We stand behind our product quality, ShelLab CO₂ Incubators come with the most extensive warranties in the industry.

Limited

- 5 years parts (labor included in the US)
- 7 years IR sensor

SC06AD SC06AD-2 High Heat Decontamination 28.8 × 30.3 × 39.3 73.1 × 76.8 × 99.7 20.2 × 20.0 × 25.5 51.4 × 50.8 × 64.7
High Heat Decontamination 28.8 x 30.3 x 39.3 73.1 x 76.8 x 99.7 20.2 x 20.0 x 25.5
28.8 x 30.3 x 39.3 73.1 x 76.8 x 99.7 20.2 x 20.0 x 25.5
73.1 x 76.8 x 99.7 20.2 x 20.0 x 25.5
20.2 x 20.0 x 25.5
51.4 x 50.8 x 64.7
3211 X 3010 X 0117
5.9
169
0-20%
+/-0.1%
<5 Minutes to 95% of Setpoint
Up to 95%
Ambient +5°C to 50°C
+/- 0.25°C at 37°C
Yes
Yes
Yes
103

*All specifications are determined by using average values on standard equipment at an ambient temperature of 25°C (77°F) and line voltages within +/-10% of unit type (115V/230V). Temperature specifications follow DIN 12880 methodology. We reserve the right to change specifications at any time.

0740519 12/13







copper shelving option

more information copper shelving

Copper Shelf Options	Part Number
Copper Shelf Slides	5820505
Copper Shelf	5820504
Copper Shelf Kit (3 Shelves, 6 Slides)	9750582



