

Accessories 2306

Rotor/Bucket

Order No.: 220.96 V04



Angle rotor 12 x 15 ml
Angle: 28°

Order No.: 220.87 V09



Angle rotor 24 x 1,5 / 2,0 ml
Angle: 45°

Order No.: 220.87 V10



Angle rotor 24 x 1,5 / 2,0 ml
Angle: 45° | hermetically sealed

Order No.: 220.97 V04



Angle rotor 6 x 50 ml
Angle: 28°

Order No.: 221.17 V03



Angle rotor 30 x 1,5 / 2,0 ml
Angle: 45° | hermetically sealed

Adapters								
Acceleration in s: 13 / 105 Deceleration in s: 13 / 300 max. Speed: 6000 rpm RCF-value: 4180 xg max. radius: 10,4 cm max. Ø: 17 mm	Volume	15 ml Falcon	15 ml	15 ml	7 ml	5 ml	1,5 ml	
								
	Tube Ø in mm	17	17	17	13	13	11	
	Length in mm	120	105-125	95-105	95-115	70-95	42	
Adapter Order No.:		without	without	701.010	701.011	701.012	701.000	

Acceleration in s: 22 / 200 Deceleration in s: 35 / 290 max. Speed: 13500 rpm RCF-value: 17319 xg max. radius: 8,5 cm max. Ø: 11 mm	Order No.:	Tube rack	Ø (mm)	Length (min./max.):
	without Adapter	1 x 2,0 ml	11	40
	without Adapter	1 x 1,5 ml	11	40
	704.005	1 x 0,5 ml	8	45
	704.004	1 x 0,4/0,2 ml PCR	6	45

Acceleration in s: 25 / 200 Deceleration in s: 35 / 340 max. Speed: 13500 rpm RCF-value: 17319 xg max. radius: 8,5 cm max. Ø: 11 mm	Order No.:	Tube rack	Ø (mm)	Length (min./max.):
	without Adapter	1 x 2,0 ml	11	40
	without Adapter	1 x 1,5 ml	11	40
	704.005	1 x 0,5 ml	8	45
	704.004	1 x 0,4/0,2 ml PCR	6	45

Acceleration in s: 13 / 105 Deceleration in s: 13 / 300 max. Speed: 6000 rpm RCF-Value: 3820 xg max. radius: 9,5 cm max. Ø: 30 mm	Volume	50 ml Falcon	15 ml Falcon	15 ml				
								
	Ø Tube in mm	30	17	17				
	Length in mm	120	120-125	95-105				
Adapter Order No.:		without	603.001	701.010				

Acceleration in s: 27 / 250 Deceleration in s: 34 / 450 max. Speed: 12000 rpm RCF-value: 15133 xg max. radius: 9,4 cm max. Ø: 11 mm	Order No.:	Tube rack	Ø (mm)	Length (min./max.):
	without Adapter	1 x 2,0 ml	11	40
	without Adapter	1 x 1,5 ml	11	40
	704.005	1 x 0,5 ml	8	45
	704.004	1 x 0,4/0,2 ml PCR	6	45

HERMLE
LABORTECHNIK

HERMLE Labortechnik GmbH
Siemensstraße 25
78564 Wehingen • Germany
Phone: + 49 7426/9622-50
Telefax: + 49 7426/9622-49
E-Mail: vertrieb@hermle-labortechnik.de
Web: www.hermle-labortechnik.de

UNIVERSAL CENTRIFUGE 2306



IVD CONFORM

Equipment

Wide Range



HERMLE
LABORTECHNIK

UNIVERSAL CENTRIFUGE 2306

Accessories 2306

The small, universal laboratory centrifuge offers a wide range of applications, which covers clinical laboratory requirements, research, as well as industrial laboratory use. The Z 306 offers many accessories, including rotors for microtitreplates and high speed fixed angle rotors for reaction vessels or centrifuge tubes up to 50 ml. Additionally, swing out rotors, which hold a volume of 4 x 100 ml can be used with special adapters for tissue culture tubes and also for all common blood tubes.

This wide range of accessories makes the Z306 laboratory centrifuge indispensable in research laboratories and small clinical hospital laboratories.




Specifications

- Maintenance-free induction drive
- Microprocessor with large LCD display
- Lid lock
- Automatic rotor identification system with over speed protection, immediately upon rotor insertion
- Active imbalance-identification and -cut off (rotor specific)
- Imbalance tolerant drive
- A large variety of rotors and accessories
- Manufactured according to international safety regulations, i.e. IEC 61010
- Audible signal at the end of each run with several melody options
- Rotor exchange is easy to handle
- Noise level under 63 dBA at max. speed
- Improved ventilation system keeps samples cool

Distinct Control Panel

- Simple one handed operation
- Easy to program with gloves on
- Splash proof digital display
- Permanent indication of Preset and Actual Values
- Selection of speed in both rpm and g-force, with increments of 10
- 10 acceleration and 10 deceleration rates, Possibility of unbraked deceleration
- Pre-selection of running time from 10 s to 99 h 59 min or continuous
- Storage of up to 99 runs incl. rotor type
- Quick-key for short runs
- Diagnostic program

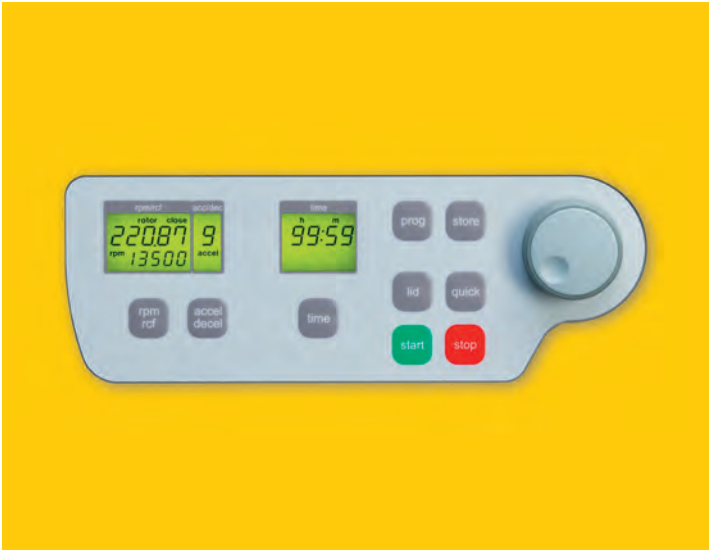

Technical Data Z 306

max. speed: 13500 rpm
max RCF: 17319 xg
max. volume: 4 x 100 ml
Speed range: 200 - 13500 rpm
running time: 59 min 50 s/ 10 s increments
99 h 59 min 50 s / 1 min increments

Dimensions (w x h x d): 35,5 cm x 33 cm x 47,4 cm
Weight: 30 kg

Order numbers:
310.00 V01 - Z 306, 230 V / 50 - 60 Hz / 240 W
310.00 V02 - Z 306, 120 V / 50 - 60 Hz / 300 W

Rotor/Bucket

Order No.: 221.12 V03



Swing out rotor 4 Places
Angle: 0-90°

Volume
Tubes Ø in mm

Length min./max. (mm)
Tubes per rotor
Radius correction (cm)
Order No.: Tube rack

Order No.: 220.72 V06

Swing out rotor 4 Places
Angle: 0-90°

605.000
Bucket for Tube Rack

605.002
Hermetic lid for bucket



Volume
Tubes Ø in mm

Length min./max. (mm)
Tubes per rotor
Radius correction (cm)
Order No.: Tube rack

Order No.: 221.16 V03



Swing out rotor 2 x 3 Microtitreplates
Angle: 0-90°

Order No.: 220.50 V08



Swing out rotor 2 x 3 Microtitreplates
Angle: 0-90°

Order No.: 221.19 V02



Angle rotor 30 x 15 ml
Angle: 35°

Adapters

Acceleration in s: 14 / 110 Deceleration in s: 19 / 250 max. Speed: 4500 rpm RCF-value: 3350 xg max. radius: 14,8 cm	Volume Tubes Ø in mm	100 ml RB Glass 44	85 ml RB/OR 38	50 ml RB Glass 34	50 ml RB Glass 29	50 ml Falcon 29
	Length min./max. (mm)	90/108	90/108	92/118	90/133	116
	Tubes per rotor	4	4	4	4	4
	Radius correction (cm)	0,7	0,7	0,3	0,3	0,3
	Order-No.: Tube Rack	626.000	626.001	626.008	626.009	626.002

50 ml Rim 29	30 ml Universal 25	25 ml RB Glass 25	15 ml Falcon 17	15 ml Rim 17	15 ml RB Glass 17	5/7 ml RB Glass 13	1,5 ml Reakt. 11	4,5 ml Sarstedt 12	5/7ml 13	5/7ml Vakut. 13
116	95/105	95/110	120	120	90/108	75/100	42	75/100	65/80	75/100
4	4	4	8	4	28	36	40	36	28	28
0,2	0,2	0,2	0,0	0,0	0,3	0,6	0,6	0,6	3,2	0,6
626.010	626.011	626.012	626.003	626.013	626.004	626.005	626.014	626.007	626.015	626.006

Acceleration in s: 14 / 110 Deceleration in s: 17 / 170 max. Speed: 4000 rpm RCF-value: 2486 xg max. radius: 13,9 cm	Volume Tubes Ø in mm	100 ml RB Glass 40	85 ml RB/OR 38	50 ml RB Glass 34	50 ml RB/OR 29	50 ml Falcon 29
	Length min./max. (mm)	90/108	90/108	92/118	90/133	116
	Tubes per rotor	4	4	4	4	4
	Radius correction (cm)	0,7	0,7	0,0	0,0	0,3
	Order-No.: Tube Rack	705.000	705.012	705.001	705.010	705.006

50 ml Rim 29	30 ml Universal 25	25 ml RB Glass 25	15 ml Falcon 17	15 ml Falcon 17	15 ml Rim 17	15 ml RB Glass 17	5/7 ml RB Glass 13	1,5 ml Reakt. 11	10 ml Sarstedt 17	5/7ml 13	10 ml Vakutainer 17	50ml Falcon 29	15ml Falcon 17
116	75/105	95/110	120	120	120	90/108	75/100	42	85/100	75/100	75/100	116	120
4	4	4	4	8	8	16	28	20	12	16	16	8	12
0,7	0,3	0,0	0,0	0,0	0,0	0,3	0,6	0,6	0,2	0,6	0,2	0,7	0,0
705.013	705.014	705.002	705.018	705.005*1	705.015*2	705.003	705.007	705.009	705.011	705.008	705.016	605.004	605.005

Acceleration in s: 25 / 230
Deceleration in s: 26 / 490
max. Speed: 4500 rpm
RCF-Value: 2740 xg
max. radius: 12,1 cm
max. Loading height: 58 mm
max. Loading weight: 310 g

Order No.: 706.000
Tube rack
48 x 1,5 ml (2,0 ml)
Tubes

Acceleration in s: 14 / 110
Deceleration in s: 19 / 330
max. Speed: 500 rpm
RCF-value: 1424 xg
max. radius: 10,4 cm
max. Loading height: 65 mm
max. Loading weight: 310 g

Order No.: 706.000
Tube rack
48 x 1,5 ml (2,0 ml)
Tubes

Acceleration in s: 19 / 157 Deceleration in s: 20 / 590 max. Speed: 4500 rpm RCF-value: 2830 xg max. radius: 12,5/10,9 cm max. Ø: 17 mm	Order No.: 701.010 701.011 701.012 701.000	Tube rack 1 x 15 ml-Glass 1 x 7 ml 1 x 5 ml (RIA) 1 x 1,5 ml	Ø (mm) 17,0 13,5 13,5 11	Length (min./max.): 95/115 95/115 70/90 42
--	---	--	--------------------------------------	--